MANUAL BREAST PUMP

with bottle 120ml

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Product Description

Easy to use and practical, LORELLI pump allows you to regulate the intensity of squeezing.

Designed for light and convenient contact with the breast, the pump is compatible both to ordinary bottles and to wide throat bottles. The pump bottle is used directly for feeding.

Product Features

- Soft shield cover fits your breast comfortably.
- Handle snugly mould into your hand.
- Easy cleaning and assembly.
- Self regulating suction power.

- The feeding teat is S-size and it is suitable only for this breast pump feeding bottle.

Sterilize:

- Prior to using your pump for the first time, disassemble it and boil all parts for 10 minutes. Remove all parts from water immediately after boiling. Completely air dry before assembly and before the next use.

- Do not wash with strong cleaning detergents. This may damage the products.

For your child's safety and health WARNING! Before and after feeding:

Before first use, clean the product.

- ALWAYS USE THIS PRODUCT WITH ADULT SUPERVISION.
- Never use feeding teats as soothers!
- Continuous and prolonged sucking of fluids will cause tooth decay.
- Always check food temperature before feeding.
- Squeeze the nipple to remove any water or milk left.
- Check the nipple before and after each feeding. Replace it with new one if it is damaged or worn out.
- Tooth decay in young children can occur even if non-sweetened fluids are used. This can occur if the baby is allowed to use the bottle for long periods throughout the day and particularly throughout the night when saliva flow is reduced.
- Heating in microwave oven may produce localised high temperatures.

Cleaning:

The manual breast pump parts may be cleaned in the upper compartment of your dishwasher or hand washed with dish soap and hot water. For best results when hand washing the pump, we recommend using the cleaning brush to thoroughly clean your breast pump. Do not use solvents or abrasives. Rinse with hot, clean water. The valve, silicone diaphragm and stem, and the Comfort cushion can be washed in the dishwasher; however, we recommend washing them by hand to prevent loss and to extend their usefulness. Completely air-dry all parts on a clean towel or drying rack.

Note: Use caution when cleaning the valve. It is important for proper suction. Do not use excessive force and do not clean with baby bottle brush, a nipple brush or the cleaning brush. Gently wash by hand in hot, soapy water and rinse under hot running water.

How to Store Your Expressed Breastmilk

Breastmilk should be stored in a clean, sealed container. Do not store your breastmilk in a freezer door. The door tends to be the warmest part of the freezer because of opening and closing.

Thawing

Thaw frozen milk in the refrigerator or place the container under cool running water, gradually warming the water until the milk is thawed. Thawed milk can be kept in a refrigerator for up to 24 hours. Do not refreeze milk once it is thawed

Gently swirl the container of thawed milk before feeding to baby to mix the layers that have separated. Breastmilk is not homogenized and the "cream" will rise to the top of the bottle. Breastmilk can vary in color and consistency depending upon the time of the day that the milk was expressed, what you have eaten, and the age of the baby at the time of pumping. Do not shake the bottle vigorously to mix the layers as doing so can damage the live components of your milk.

Do not microwave human milk! Microwaving alters the composition of the milk. Do not thaw breastmilk in very hot or boiling water.

Instructions for safe feeding

- Always wash your hands before preparing the food and before feeding the baby.
- Before using the product always wash and sterilize.
- Do not leave the nipple on direct sunlight, near to heat sources or in sterilizing solution for longer than recommended.
- Keep all parts of the product out of children's reach!
- Never dip the nipple in sweet substances! This may cause tooth decay. Replace the nipple each one or two months to keep its properties and safety!
- Always wash after each feeding.

Recommendations:

Drawing milk from breasts can be done for many reasons:

- If your breasts produce too little milk and you want to increase the production: using breast pump stimulates feeding process, which leads to lactation increase.
- After breast feeding when you still have milk.
- When you notice sudden milk overfilling in your breasts (milk lots), pumping it out will loosen them and give you a relief.
- If your nipple goes inside the overfilled breast drawing a little milk before the actual feeding time will soften the areola of the mamma (edging) and bring the nipple forth, by which it will make it easier to grab.

When you feel sudden overfilling between feeding times, you can deal with it without waking up the baby.

No Milk Expressed: If you are not successful in expressing any milk, try the following:

- Make sure that all the parts have been properly dried and assembled and that proper suction is created.
- Try to relax as successful pumping is a learned art and may take several practice attempts to get started.

Pain When Pumping - If you are experiencing pain when expressing milk, try the following:

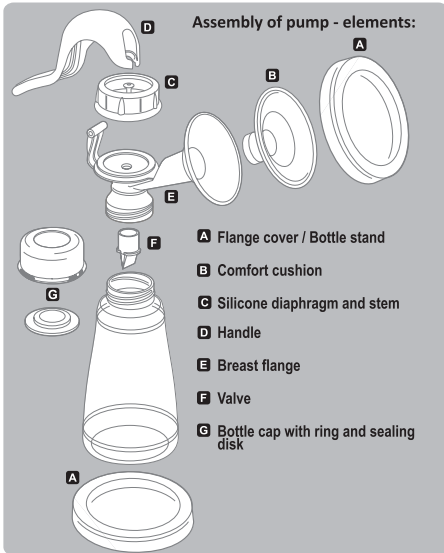
You may be pumping too hard. It may not be necessary to use all of the suction that the pump can generate. Try depressing the handle half way and remember the 2-3 second suggested rhythm.

Breastmilk Storage Guidelines

Where	Temperature	Time
At room temperature (fresh milk)	66 ° to 78 ° F (19 ° to 26 ° C)	4 hours (ideal) up to 6 hours (acceptable)* (Some sources use 8 hours)
In a refrigerator	<39°F (>4°C)	72 hours (ideal) up to 8 days (acceptable)**
In a freezer	-0.4° to -4°F (-18° to -20°C)	6 months (ideal) up to 12 months (acceptable)

* The preference is to refrigerate or chill milk right after it is expressed.

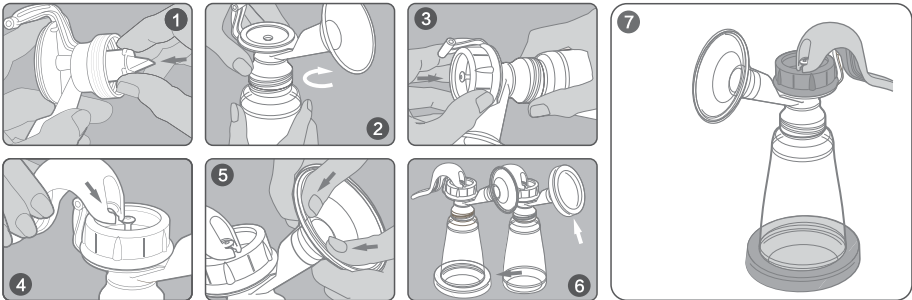
** Eight days is acceptable, ideally collected in a very clean, careful way.

**How to assemble your breast pump:**

Note: Prior to assembling your Breast Pump ensure that all parts have been cleaned and sterilized. Pictures on the page show a detailed instruction of how to dismantle and put back together the breast pump. Wash your hands thoroughly.

1. Gently secure the valve onto the breast flange from underneath, making sure the valve is flush with the lip of the breast flange. It is not necessary to use excess force when applying the valve. If you find that the valve is difficult to remove, you have pushed it too far onto the breast flange.
2. Screw the bottle onto the breast flange.
3. Place the silicone diaphragm and stem on top of the breast flange. Make sure it fits securely around the rim by pressing down with your fingers to ensure a perfect seal.*
4. Place the forked end of the handle under the diaphragm stem and push down gently on the handle until it clicks into place.
5. Carefully insert the Comfort cushion into the funnel of the breast flange, making sure it is perfectly sealed around the rim of the funnel.*
6. To make the pump completely stable, insert the bottle into the bottle stand. Or, when not in use, the bottle stand can be used as a cover over the breast flange (flange cover).
7. When fully assembled, this is what your pump should look like.

** These parts are designed to fit snugly. For easier assembly, just dampen the rim of the silicone diaphragm and/or Comfort Cushion with water.*

**How to use your breast pump:**

Before starting to use the pump, wait until your milk supply and breastfeeding schedule are well established. This is normally at least 3-4 weeks after birth.

1. Prepare clean, sterilized and ready to use breast pump.
2. Wash hands with soap carefully.
3. Warm the breast before the process: bath towel or warm shower.
4. Massage breast with circular moves and slightly rub the nipple, just like during breast self-examination. This stimulates milk production as well as its flow.
5. Place the funnel of the pump on the breast so that the nipple gets inside it and no air is coming through. There must be vacuum for the pump to work properly. If the inside of the funnel gets damp, the nipple will freely translocate.
6. It is advised to drink a glass of water during the procedure.
7. Milk obtained after every session should be placed in separate bottles.

WARNING! Should you have any doubts whether you store your milk properly, it will be safer for your baby if you contact the doctor or feed the baby with freshly drawn milk.

8. As you gently push down on the pump handle you will feel the suction on your breast. You do not need to depress the handle fully to create a vacuum; only as much as is comfortable. Initially squeeze the handle with short quick compressions until let-down is achieved. Your milk will soon start flowing even though you may not use all the suction the pump can generate.
9. Begin by pumping 6-7 times quickly to initiate let down. Then, hold the handle down for 2-3 seconds, and allow it to return to its resting place. These 2-3 second cycles imitate your baby's natural nursing pattern and allow the milk to flow into the bottles between strokes.

10. Milk should begin to flow after the first few strokes. Do not worry if your milk does not flow immediately. Relax and continue pumping.

Warning! Do not continue pumping for more than five minutes at a time if no result is achieved. Try expressing at another time during the day. Should pumping become very painful or uncomfortable, you should stop using the pump and consult your doctor.

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Producer: Didis Ltd.